anti-ZKSCAN8-RAB-C51



Applications

Competition ELISA	Western Blot	SPR	Spiked IP	Immunofluorescence	IP-MS	ChIP
Pass			Pass	Pass		

*rAb has been tested for the following applications. See below for the experimental details.

Antibody information

rAb ID: anti-ZKSCAN8-RAB-C51


Description: recombinant Fab fragment obtained by recombinant antibody (rAb) phage display recognizing *ZKSCAN8* protein under non-denaturing conditions; specificity and affinity tested.

Binder type: rAb **Isotype:** IgG1 **Species:** *Homo sapiens* **Produced in:** *E. coli* **rAb tags:** no tag; no tag **Specificity:** reacts with *Homo sapiens* ZKSCAN8 **Epitope:** binds to folded domain amino acids 45-132 **Storage conditions:** short term – store at n 4°C (over 6 months), long term - PBS -20°C or -80°C **Link:** http://recombinant-antibodies.org/binders/anti-ZKSCAN8-RAB-C51

Antigen information


Protein Name: Zinc finger protein with KRAB and SCAN domains 8 HGNC Symbol: ZKSCAN8 HGNC ID: 12983 Species: Homo sapiens UniProt AC: Q15776 UniProt KB: ZKSC8_HUMAN Protein Sequence: MSGLNDIFEAQKIEWHEHHHHHHENLYFQSHMLGQEVFRLRFRQLRYQETLGPREALIQLRALCHQWLRPDLNTKEQILEL LVLEQFLTILPEELQTLVKEHRLENGEEEVTLLEDLERQ Tag N-terminus: MSGLNDIFEAQKIEWHEHHHHHHENLYFQSHM Tag C-terminus: Vector Type: pET15Avi6HT_NESG Vector Link: http://beta.labgeni.us/registries/DNASU/pET15Avi6HT_NESG/ Protein Sequence Position: 45-132 Antigen source: *E. coli* Source Lab: Rutgers Source Lab ID: HR8500A.002 Description: affinity purified recombinant protein

Validation data

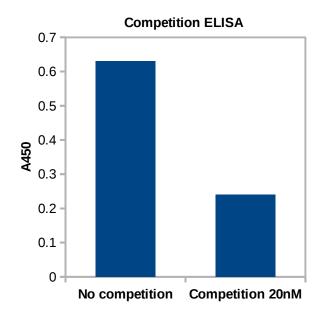

Single point competition phage ELISA Plot represents specific binding of a target to the rAb-phage in solution (right bar) in comparison to binding to the target immobilized on the plate surface (left bar). Experimental conditions were calibrated to capture binders with dissociation Constant (K_D): 20nM or lower.

Experimental Conditions: Culture supernatants containing rAbphage were diluted five-fold in phosphate-buffered saline, 0.5% (w/v) BSA, 0.1% (v/v) Tween 20 either with or without soluble antigen competitor at 20 nM. After 1 h incubation at room temperature, the mixtures were transferred to neutravidin coated plates preloaded with 50 µL of 20 nM biotinylated antigen and incubated for 15 min. The plates were washed with phosphatebuffered saline, 0.05% (v/v) Tween 20 and incubated for 30 min with horse radish peroxidase/anti-M13 antibody conjugate (1:5000 dilution). The plates were washed, developed with 3,3',5,5'-Tetramethyl-benzidine/H₂O₂ peroxidase substrate (Thermo Scientific), guenched with $1M H_3PO_4$, and the absorbance at 450 nm (A450) was determined.

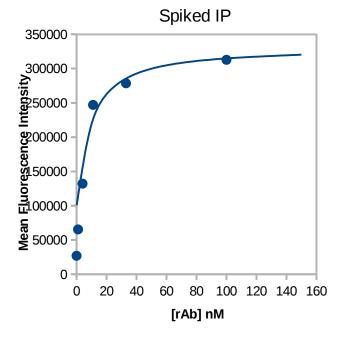

Spiked IP Tritration curve of rAb against antigen of interest. The K_D values were obtained by the least-squares fitting of fluorescence saturation data.

Dissociation Constant (K_D): 6 ± 1 nM

Experimental Conditions:

Spiked IP: Antigen was immobilized to M280 Dynabeads. A rAb, 50 nM, was pulled down from high salt AFC buffer with or without HEK293 lysate (OD280 ~10). Beads were washed with low salt AFC buffer, and the captured antibody was quantified with an anti-Fab fluorophore labeled antibody on a flow cytometer.

Affinity Measurement: Antigen was immobilized to M280 Dynabeads and incubated with a rAb, in varying concentration (100 nM down to 1 nM in three-fold dilutions). Beads were washed with BSET/BSA and quantified.

Buffers:

High salt AFC buffer: 10 mM Tris-HCl, pH 7.9, 420 mM NaCl, 0.1% NP-40

Low salt AFC buffer: 10 mM Tris-HCl, pH 7.9, 100 mM NaCl, 0.1% NP-40

PBSE/BSA: 20 mM Na2HPO4, pH 7.5, 150 mM NaCl, 1 mM EDTA, 0.5% BSA

PBSET/BSA: PBSE/BSA + 0.1% Tween-20

Immunofluorescence:

Status: Pass

Experimental Conditions: http://recombinant-antibodies.org/protocols/immunofluorescence

IP-MS – immunoprecipitation for mass spectrometric analysis:

Status:

Experimental Conditions: http://recombinant-antibodies.org/protocols/IP-MS

ChIP – chromatin immunoprecipitation:

Status:

Experimental Conditions: Pending

IP – immunoprecipitation:

Status:

Experimental Conditions: Pending

SP Elisa:

Status:

Experimental Conditions: http://recombinant-antibodies.org/protocols/ELISA-IC50-EC50-direct-coating

Visit us at http://recombinant-antibodies.org/

Contact: Recombinant Antibody Network admin@recombinant-antibodies.org

The University of Chicago

Knapp Center for Biomedical Discovery Rm. 3240G 900 E. 57th St., Chicago, IL 60637 Phone: +1 (773) 834-2776

University of California, San Franciso Byers Hall Rm. 503 1700 4th St., San Francisco, CA 94158 Phone: +1 (530) 341-2371

University of Toronto

Best Institute Rm. 117 112 College Avenue, Toronto, Ontario M5G 1L6 Phone: +1 (416) 978-1594